



NIPRO MEDICAL INDIA PVT. LTD

CUSTOMER COPY

Mohali branch Address:

NIPRO MEDICAL INDIA PVT.LTD, khatan No 761/786, Chandigarh
Zirakpur Road, Village Bhaghat Zirakpur, Mohali, Punjab Phone :- 18002095200

INSTALLATION/SERVICE REPORT

Report No: **050**

Hospital Name : *Civil Hospital*
Address : *Jagdhari (Desa)*
Nephrologist:
Chief Technician:
Biomedical Engineer:

Model : *55 Plus*
Serial No : *18CX10426, 27, 28, 29, 31, 32, 18CX1197*
Nature of Call : *pmg*
AMC/CAMC/Warranty/On call : *18CX12442*
Installation date : *09/02/2019* *19CX14576*
Warranty Expiry date : *08/02/2024* *77*
Hour Meter Reading : *9708, 9790, 3026, 6996, 9792, 5312, 6405, 7509, 856*

Errorcode: *Preventive Maintenance Service*

Machine Software Ver : Type :

Action Taken *Done Preventive Maintenance Service of II*
No.s of Serial 55 Plus M/C's.
Also done dialysate Calibration. Temp. → 37.0C
Concl. - 314.0ms/cm


Parts to be Used/Replaced: Material code Description: Qty :

Call status completed Yes No

Comments (Service Engineer): *Now all M/C function properly*

Comments (Unit-In-charge):

Signature of Service Engineer *Ravinder K...* Seal & Signature of Customer *[Signature]*
Date: *08/05/22* Time: Date: *8/5/22* Time:

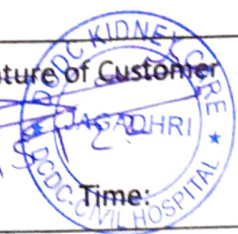

- Please follow standard disinfection protocol as recommended.
- AAMI standard RO water should be provided to Hemodialysis machines.
- It's mandatory to run the machine only with servo stabilizer or UPS.
- If spare part damaged due to power fluctuation while not using UPS or servo stabilizer are not covered under